Supporting information:

Mechanism of Ion Transport by Fluorescent Oligoester Channels

Joanne M. Moszynski and Thomas M. Fyles*

Department of Chemistry, University of Victoria. Victoria, BC. Canada. V8W 3P6

tmf@uvic.ca

Summary bilayer conductance for HO₂C-Hex-Dip-Hex-Hex-OH

Bilayer clamp experimental details.

The methodology was derived from that previously described (1). A model BC-525A bilayer clamp (Warner Instrument Corp.) was used for all planar bilayer experiments, ClampEx 8 and ClampFit 10 (Axon Instruments) were the software used for acquisition and analysis, respectively. Cups used were made of polystyrene and had 250 μ m diameter apertures. The lipid used in all cases was diphytanoylphosphatidylcholine (diPhyPC) (Avanti Polar lipids). A stock solution of 25mg/mL lipid in CHCl3 was dried under N2 and then re-suspended in 200 μ L decane. For compounds that had to be preloaded into the lipid, 0.1 – 1mol% compound in CHCl3 was added to the lipid mix and then dried down. The electrolytes used were 1M KCl, NaCl, CsCl, or NMe4Cl in 10mM HEPES, 10mM TRIS, pH 7, (unadjusted).


The aperture was primed with 0.5-1µL of decane/lipid, excess solvent was removed by blowing N2 over the aperture. The cup was then placed into the electrolyte-filled holding cell, consisting of 5mL and 3mL chambers, and salt bridges (KNO3/Agar) and electrodes (Ag/AgCl) were attached. Bilayers were formed by brushing on 1- 1.5µL of the decane/lipid mix over the aperture, and were monitored for stability, capacitance and resistance for at least 20 minutes before test compound was added. Test compounds were added by injection from an organic solution (typically no more than 1-10µL of solution) or by breaking the lipid-only bilayer and brushing on the compound-preloaded in a lipid mix. All data were hardware filtered (8- pole Bessel filter, 1 kHz) and data was collected in a survey mode using the Gapfree protocol. Bilayers were tested repeatedly for capacitance and resistance. Once formed, 'activity' from pristine bilayers was never observed.

Conductance data summary format.


Each experiment conducted is summarized in a separate file within this portfolio. Use the Bookmarks for navigation as the pages are not numbered. The structure of each summary is similar consisting of:


- A title of the form: file identifier {date} [compound] details of the experiment
- A plot of the applied potential as a function of the overall experiment time. Color overlays on this plot show the durations of the recorded traces, numbered sequentially 0000 00mn. The plot also shows when bilayer breakage occurred (if it did).
- A bar graphic showing the composition of the experiment as it proceeded. The lines establish that the indicated component was continuously present.
- Activity grid summaries for the individual traces recorded. These grids use the notation presented in ref (2). The horizontal position gives event duration in log (duration) from 1 to 10⁵millieconds. The vertical position gives conductance in log (conductance) from 1 to 3000 pS in 0.5 log unit steps (1,3,10,30etcpS). The colour codes are: green = regular "square-top behaviour; yellow = "flicker" behaviour; blue = "multiple open states" behaviour; red = "spike"

- behaviour; purple = "erratic" behaviour; grey indicates the region cannot be observed due to the instrument and data processing settings of the experiment .
- For each trace recorded there are one or more pages showing the overall trace and selected sections expanded to highlight the various conductance behaviours noted on the activity grid for the experiment.
- (1) Eggers, P.K.; Fyles, T.M.; Mitchell, K.D.D.; Sutherland, T., J. Org. Chem. 2003,68, 1050-1058.
- (2) Chui, J.K.W.; Fyles, T.M., Chem. Soc. Rev. 201241, 148-175.

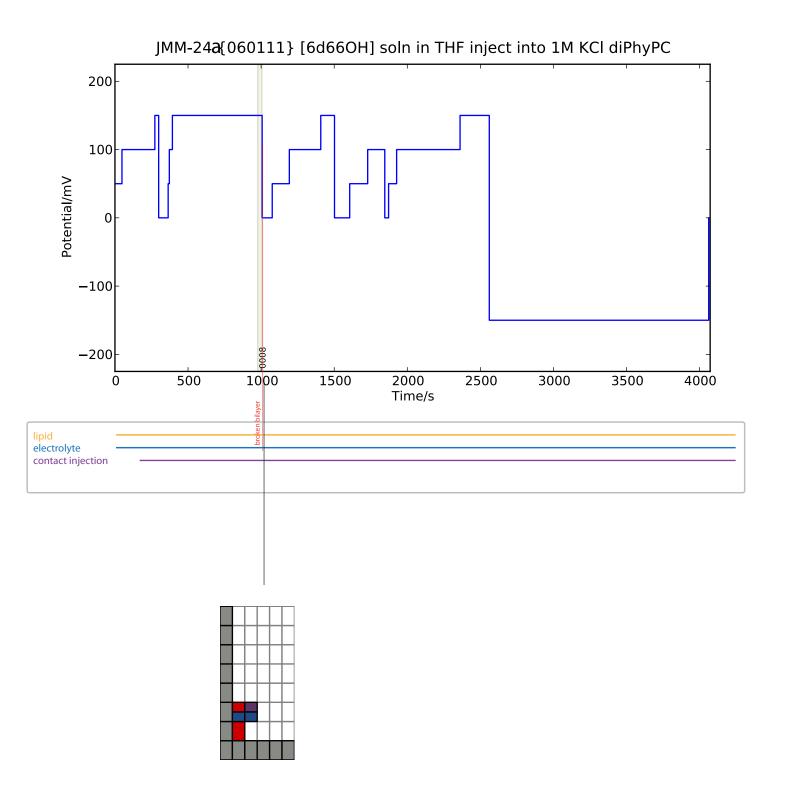



JMM-24a [6d66OH] KCl 0008

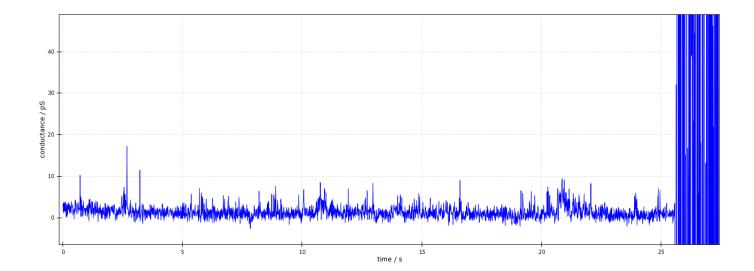


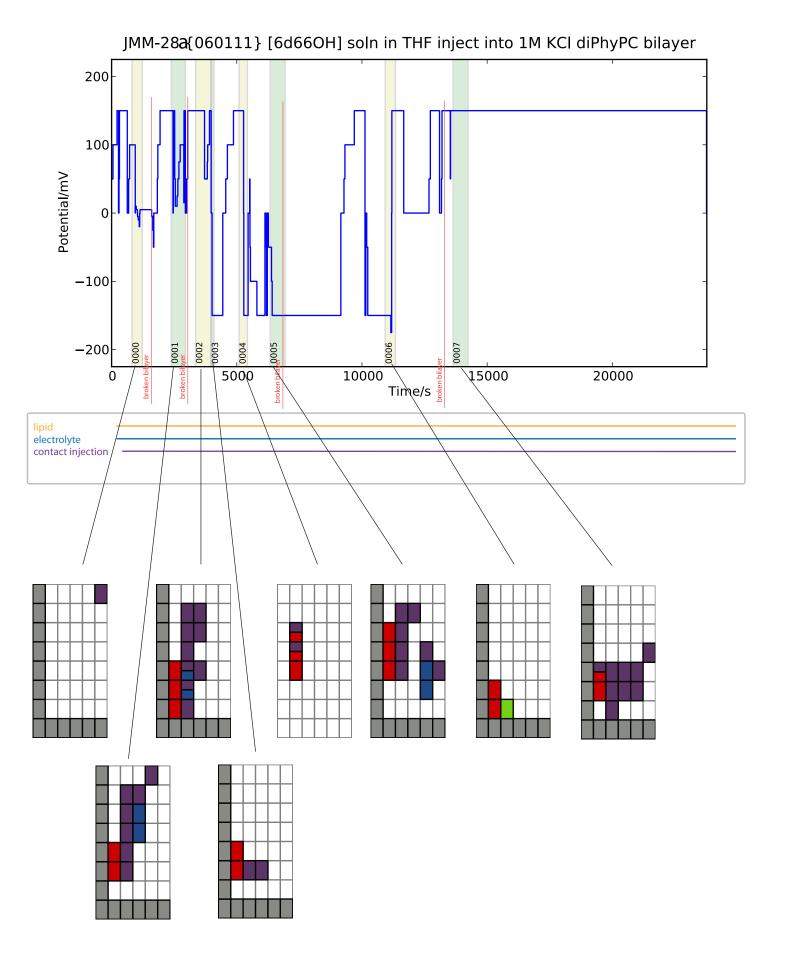



JMM-28a [6d66OH] KCI 0000

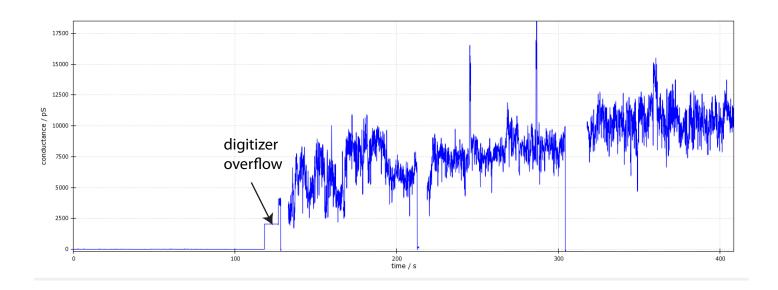


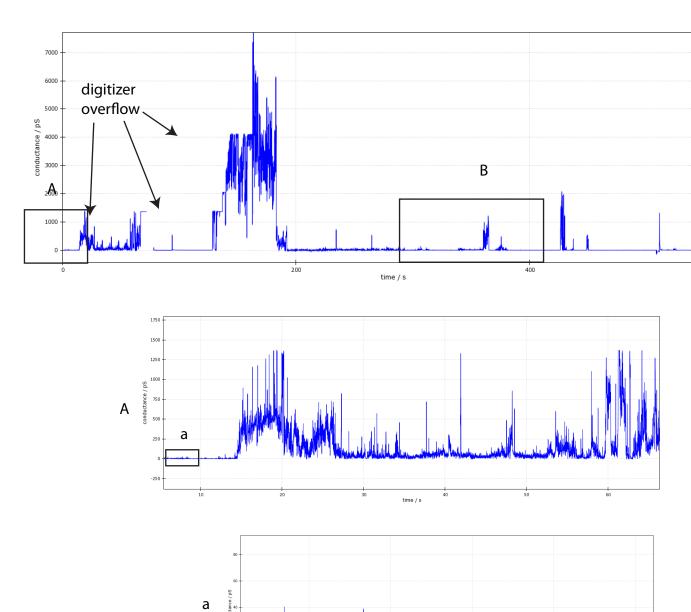

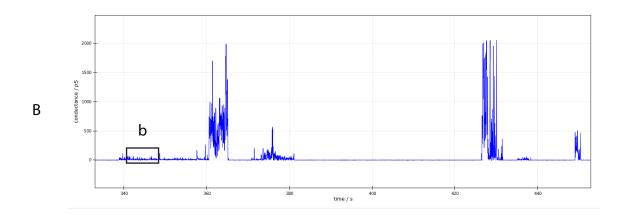


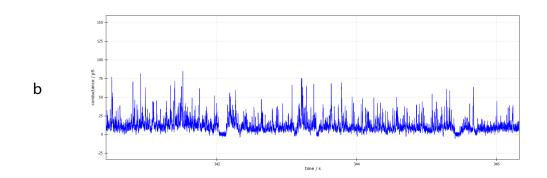


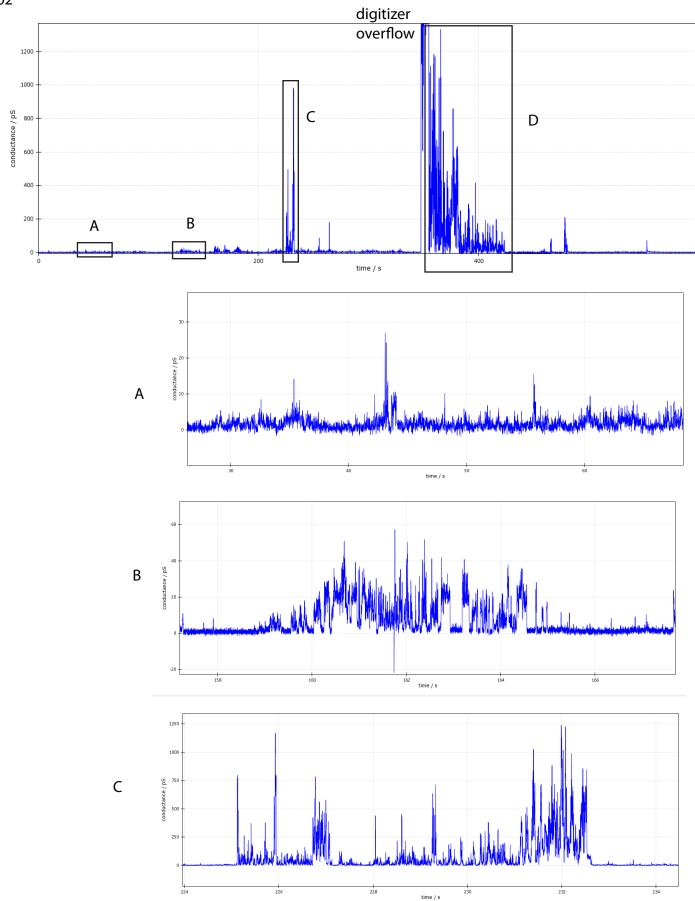


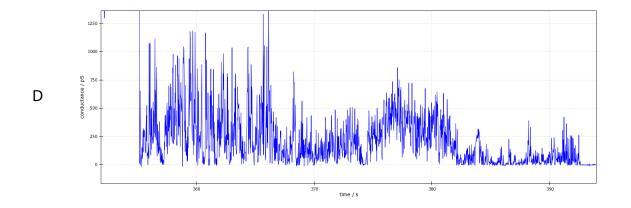


JMM-28a [6d66OH] KCI 0003

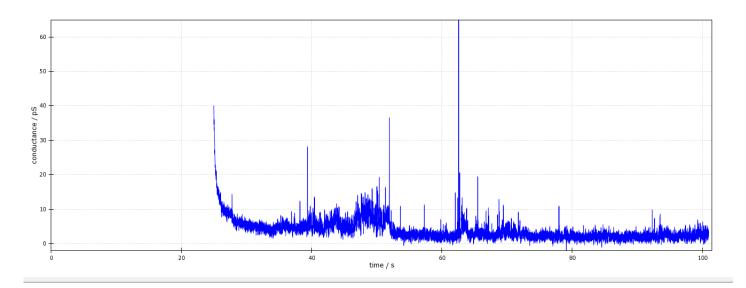


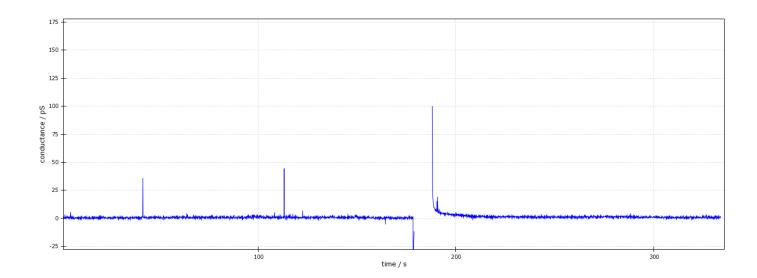


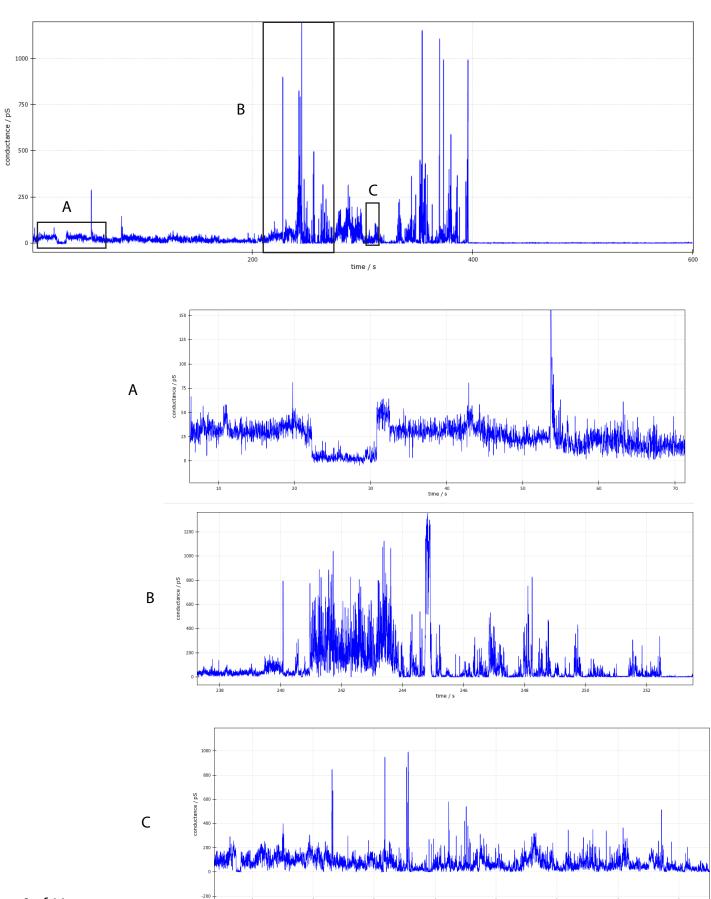



time / s

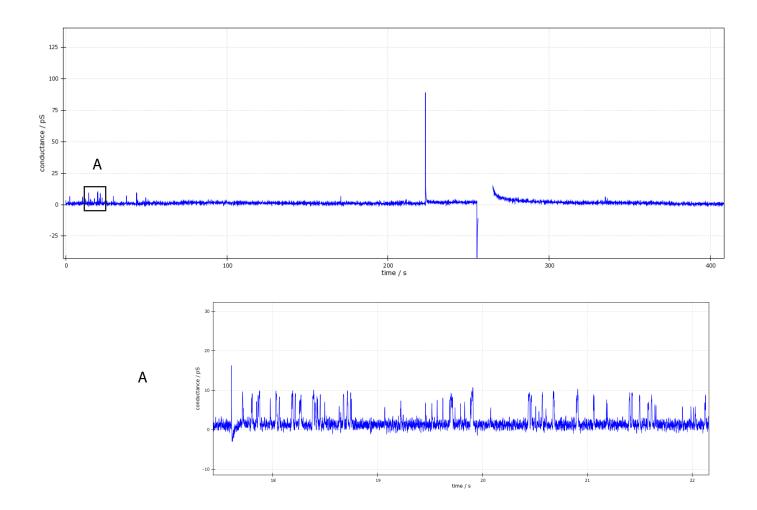


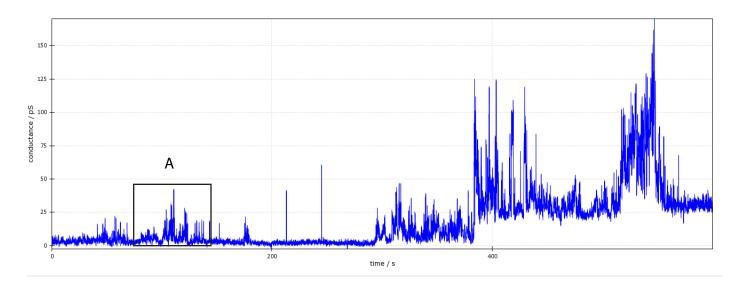


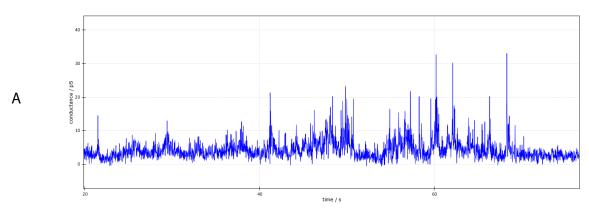


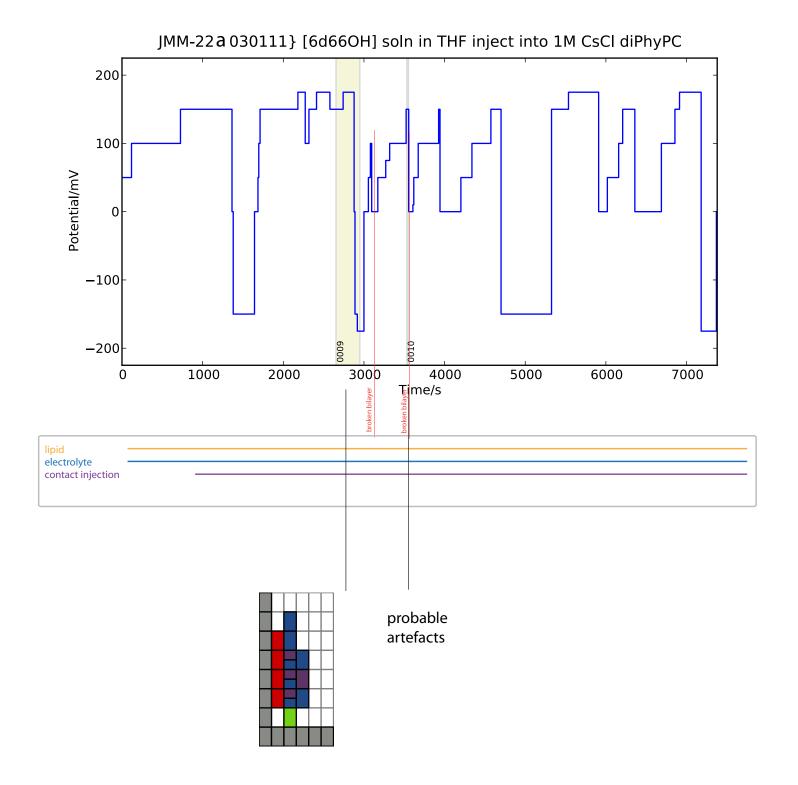


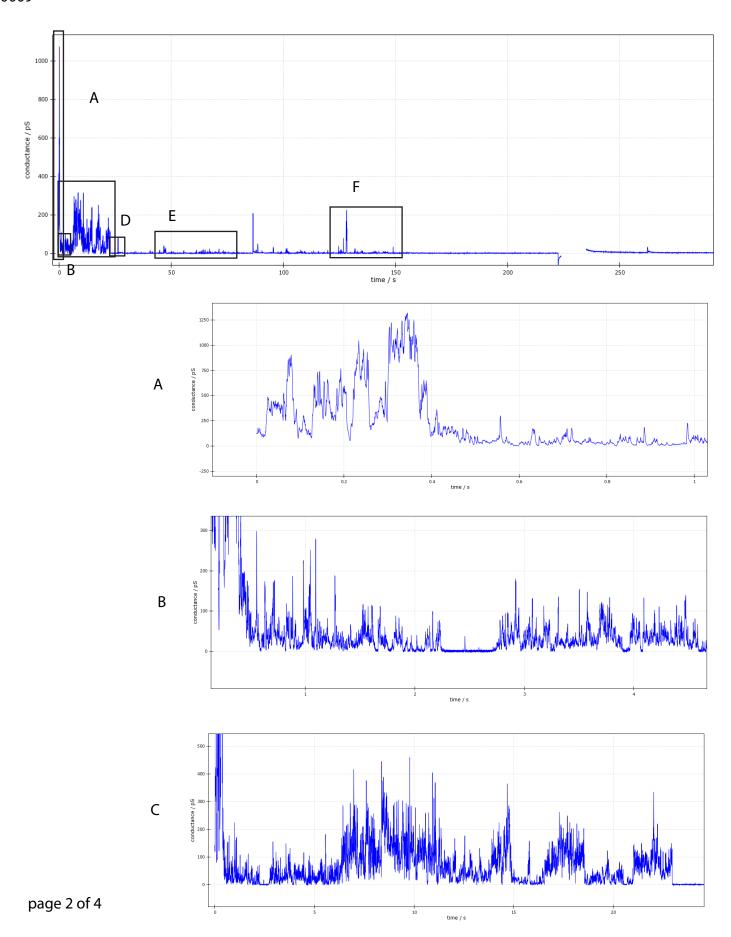


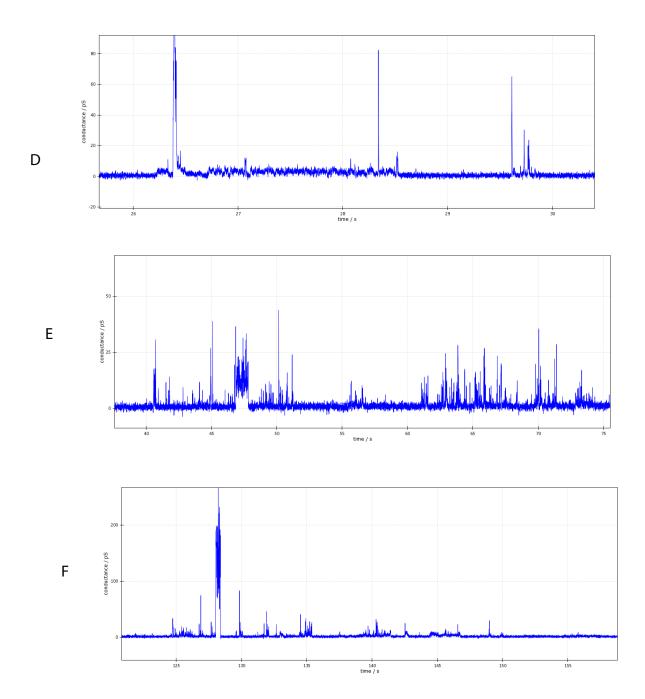


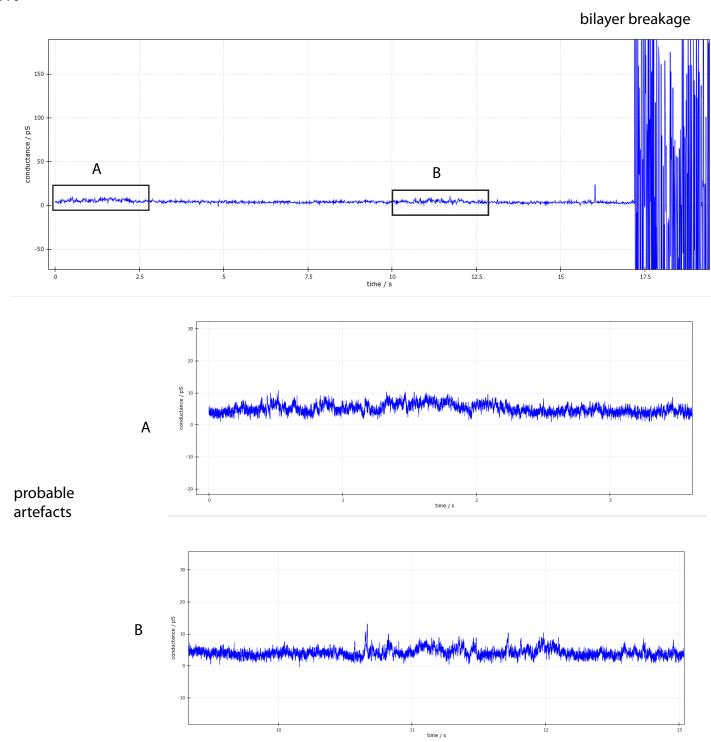


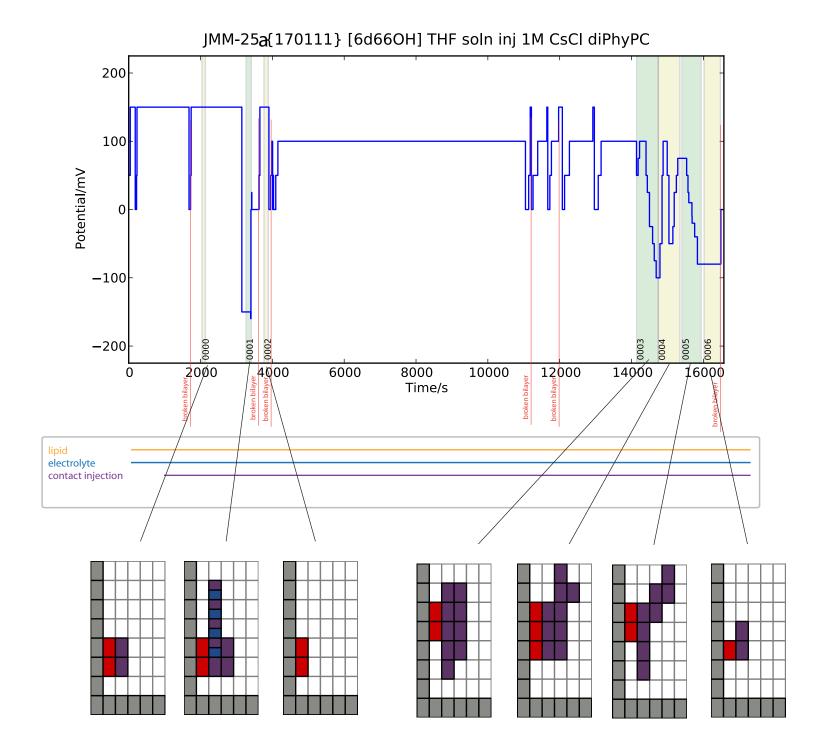


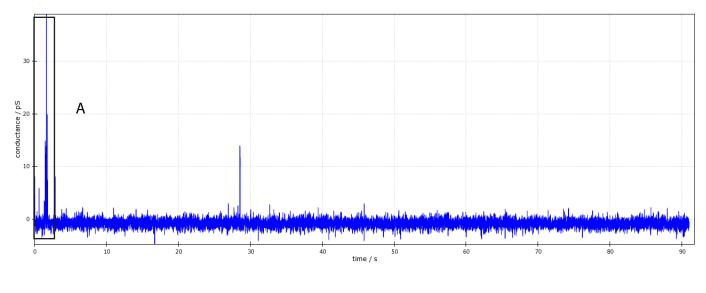


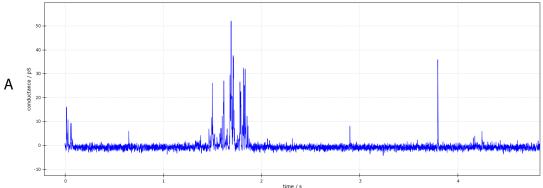


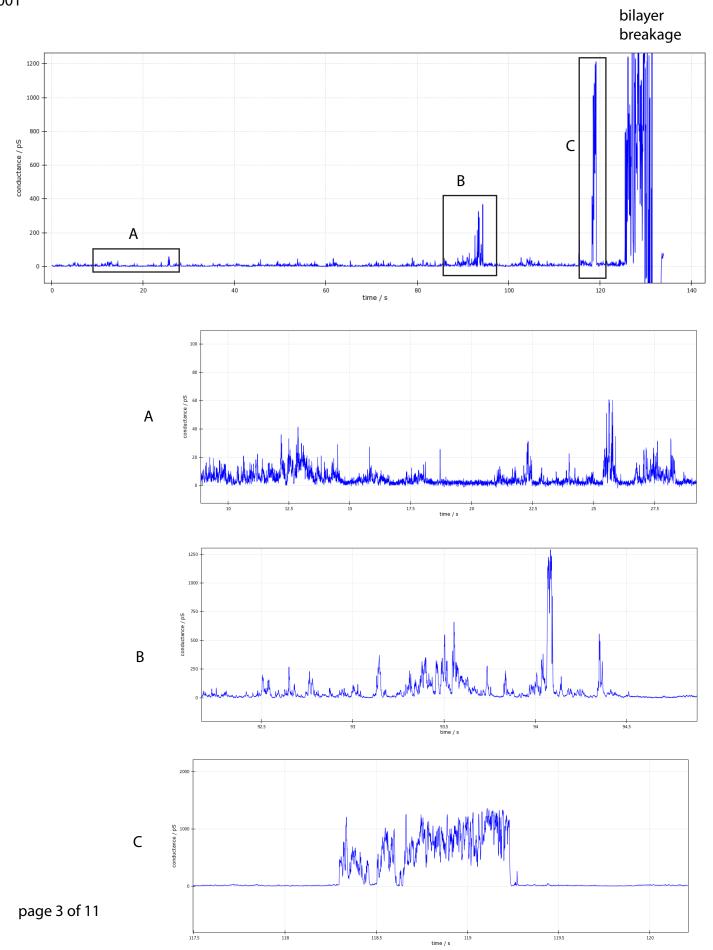


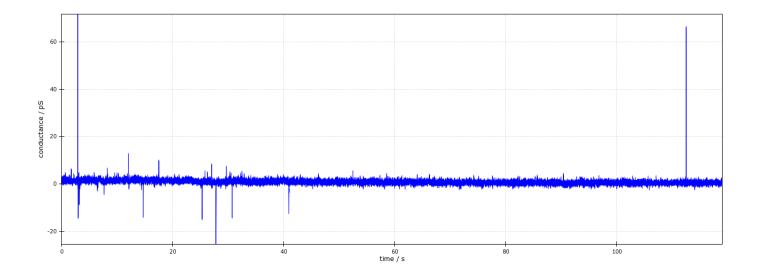


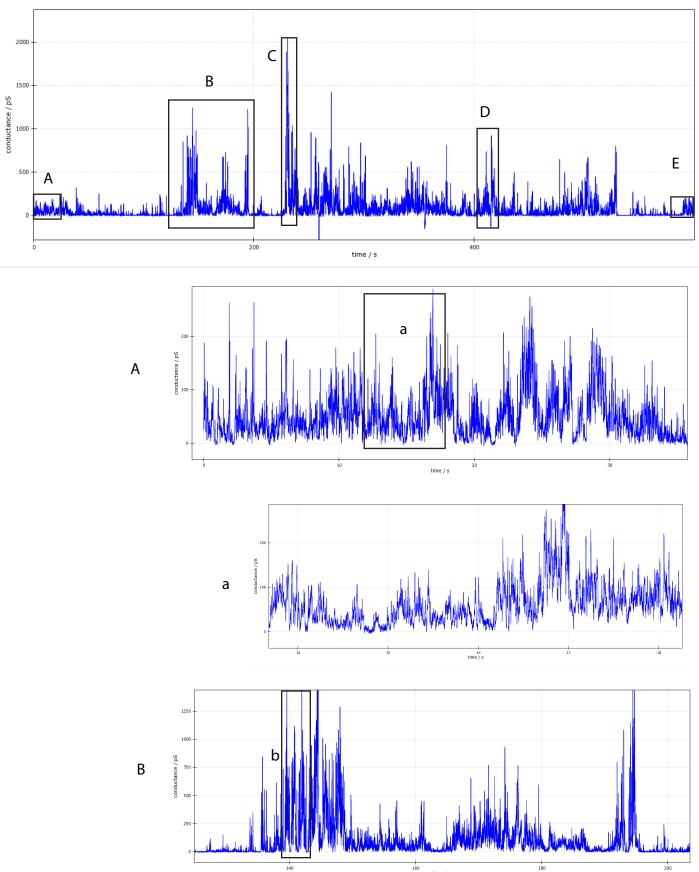

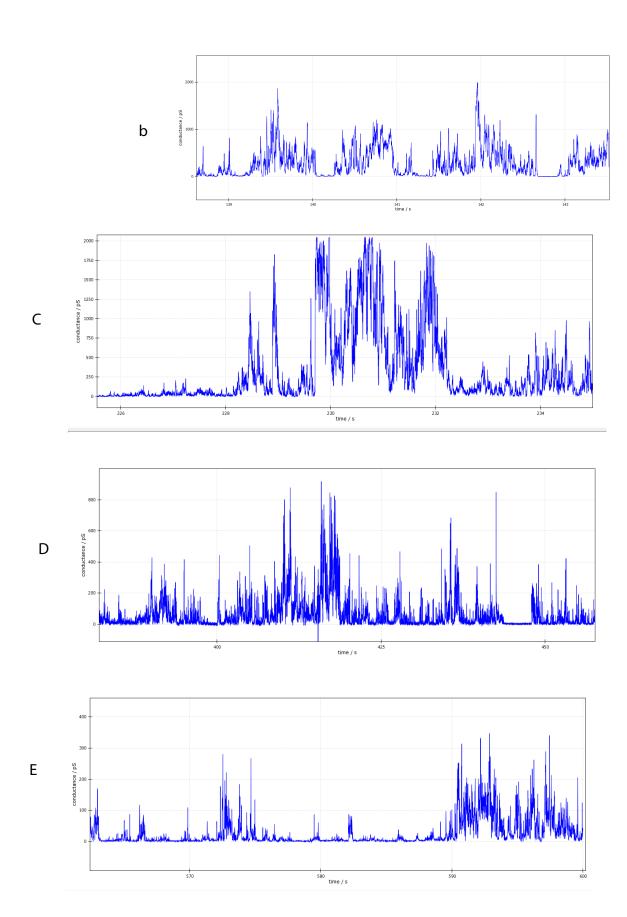

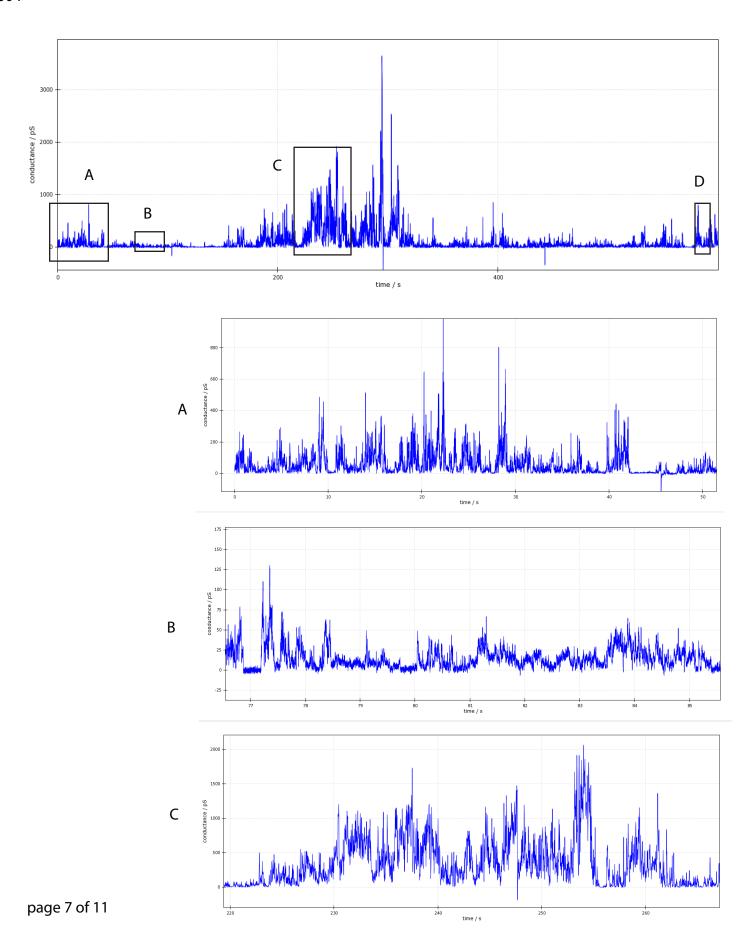

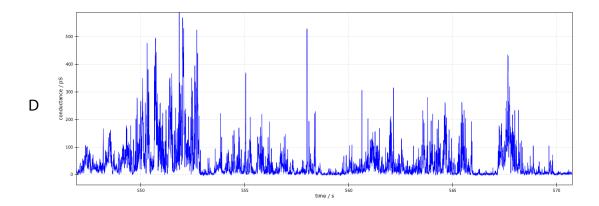

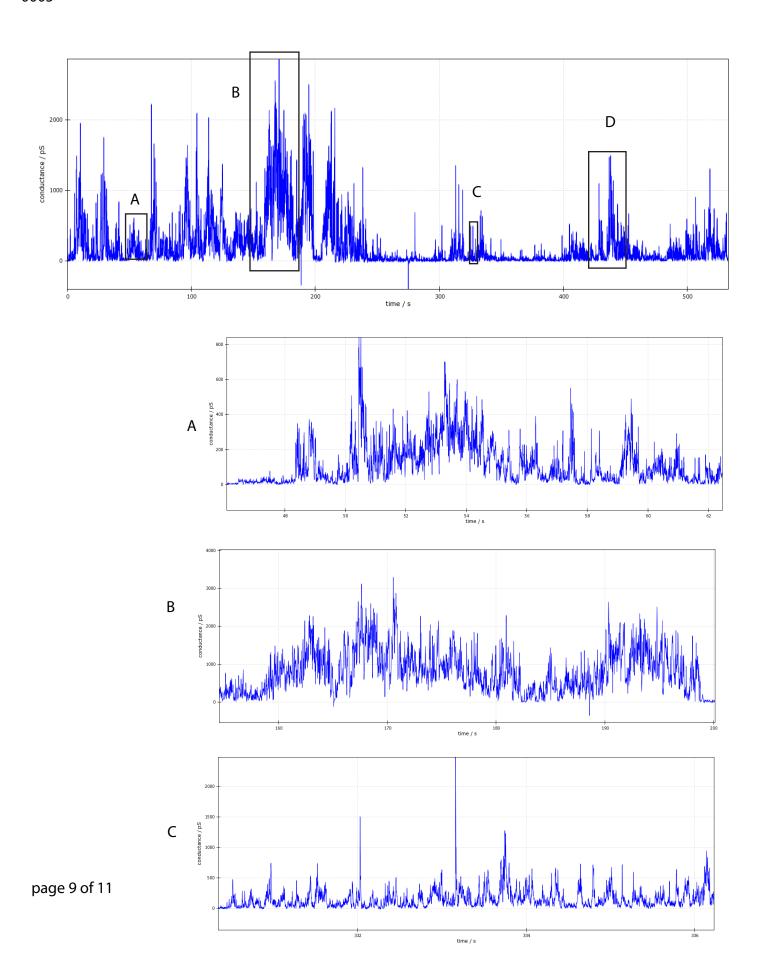

JMM-25a [6d66OH] CsCl 0005 continues

